

FLUORO-TRACK STRETCHER



FLUOROSCOPY AND C-ARM ACCESS STRETCHER

493RPCST FLUORO-TRACK STRETCHER


Designed to offer the security of a treatment table with the mobility of a stretcher. The Fluoro-Track Stretcher is designed for use with C- Arms for fluoroscopy-guided procedures. Its unique design, with a proprietary sliding top, enables a broader range of imaging capacity without repositioning the patient.

SLIDING TOP:


- Easy to operate sliding top enables a broad range of imaging capacity without repositioning the patient, enabling shorter procedure times
- Full length, radiolucent tops translate to either end providing excellent patient access
- 18 inch total range top slides 12 inches (30.5 cm) from center position toward foot end and 6 inches (15.2 cm) toward head end



slides 6 inches (19 cm)

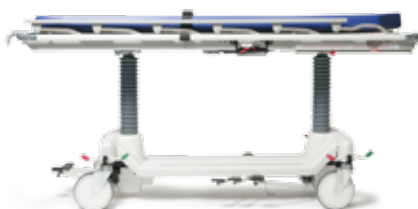

centered

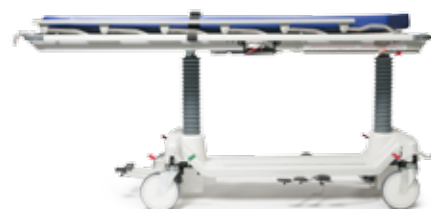

slides 12 inches (30 cm)

FEATURES:

- Low profile trimmed base creates an excellent window for C-arm access from both sides
- Designed to be a procedure table with the mobility of a stretcher, it provides improved patient load ratio for easier transport reducing risk of patient and staff injury with patient transfers
- Dual side foot pedal hydraulic controls for easy height adjustment
- Fold away side arm rails with built-in push-steer handles
- Dual column base provides superior support for the patient
- 8 inch Tente® casters facilitate maneuverability
- 3 inch comfort pads

★ See Chart for Options and Specifications



Specifications:		493RPCST	
		inch	cm†
Overall Width:		29.75	76
Overall Length:		80.25	204
Overall Height: <i>range - minimum</i>		23	58.4
Overall Height: <i>range - maximum</i>		36	91
<i>patient surface</i> Width:		24	61
<i>patient surface</i> Length:		76	193
<i>upright position</i> Retracto™ Rail Length:		64	162.5
<i>upright position</i> Retracto™ Rail Height:		15	38
<i>mid-position</i> Retracto™ Rail Height:		12.75	32.4
C-Arm Clearance		23	58.5
<i>full length</i> Radiolucent Pad:		3	7.6
<i>easy rolling</i> Tente® Casters:		8	20.3
		degree of movement™	
Backrest: <i>movement</i>		0 - 80	
<i>standard</i> Trendelenburg:		0 - 15	
<i>reverse</i> Trendelenburg:		0 - 15	
Weight Capacity: EVENLY DISTRIBUTED	lb	kg†	
	500	227	
Shipping Weight:	364	165	
Operation:		Hydraulic	

■ Standard
 ■ Optional
 — Does Not Apply
 † standard rounding method used to show metric values as whole numbers

[493 Product Info](#)

Color variations may occur due to monitor or printer settings.



© 2022, GF Health Products, Inc. All Rights Reserved.
Hausted is a registered trademark of GF Health Products, Inc.



+1 770.368.4700

Information contained herein is subject to change without notice.
The most current and complete product information can be found on our website.
www.hausted.com | www.grahamfield.com

graham field

GF1400052 Rev C22

Product Hyperlinks Enabled