


# Silkann<sup>®</sup>

# PERFECTION

for injections

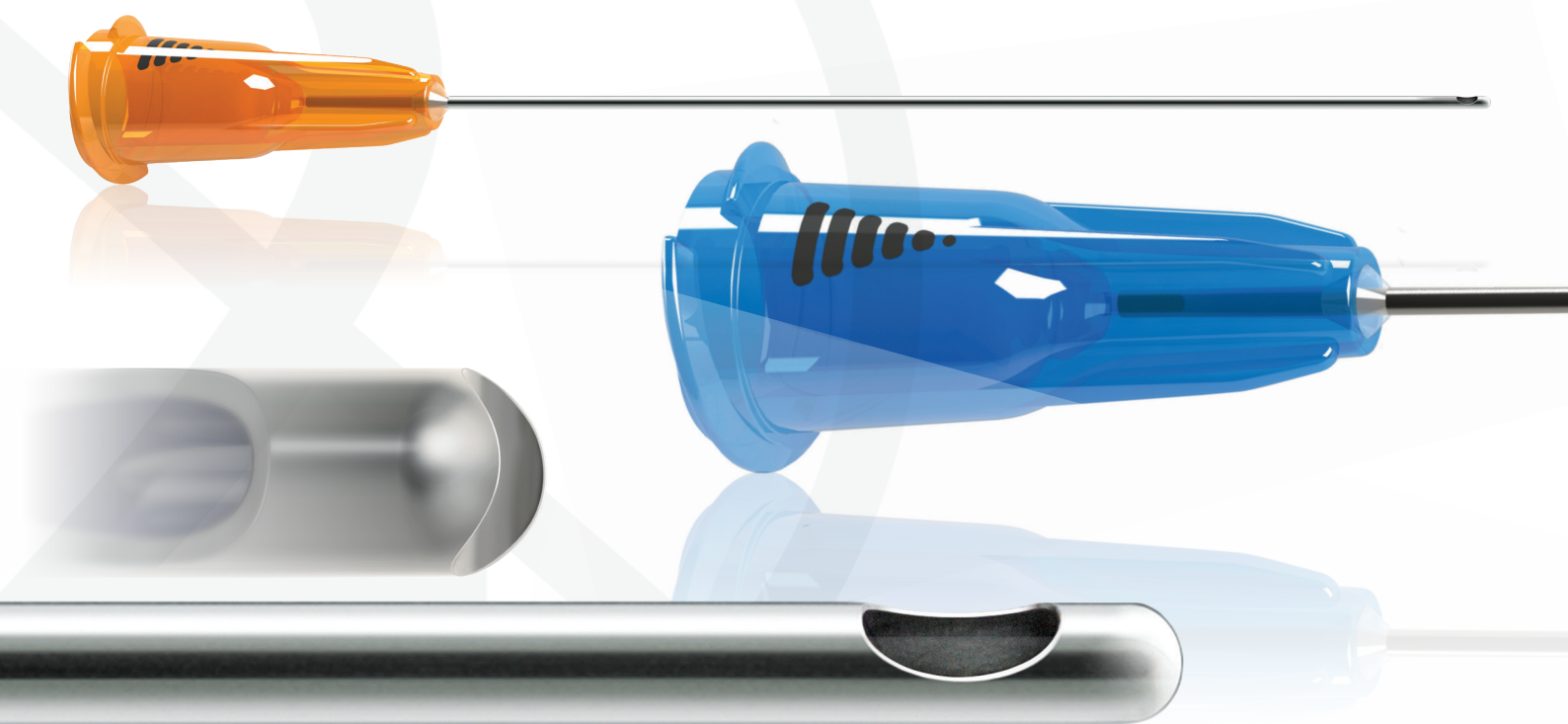

**Sterimedix<sup>®</sup>**  
aesthetic

# Silkann<sup>®</sup>


Aesthetic injections  
need a cannula that is reliable,  
consistent and cost effective:  
**perfect in every way.**

Silkann cannulas are made to the highest standards to guarantee a user experience that assists in the best possible outcome for your patients.

Design is everything. Silkann cannulas are made to ensure that injection forces are ideal for the user and patient. A steady hand allows for 100% focus on the procedure.



Sterimedix ultra thin cannula walls increase the inner diameter of the cannula and reduces injection forces

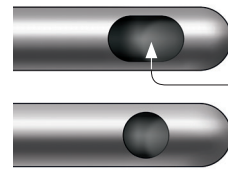


Sterimedix port hole shape matches the cannula diameter to optimise control and filler flow rate.



**Silkann cannulas feature a unique polycarbonate hub, with a screw thread to reduce the risk of detachment under pressure which is compatible with all types of filler syringes.**




- Indicator on the hub to show the port orientation
- High quality packaging to ensure the integrity of the products
- Manufactured in the UK by Sterimedix
- Compliant with all national and international standards for medical devices.

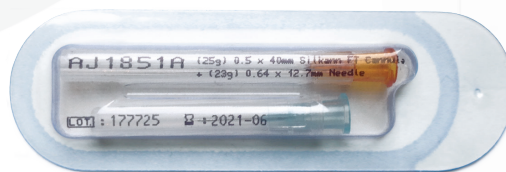All Silkann products are conveniently packed in a blister containing both the cannula and pre-hole needle





## DERMAL FILLER CANNULAS

### 22 GAUGE CANNULA + 21 GAUGE PRE-HOLE NEEDLE


			
AJ 1831A	22g TW x 70mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 21g pre-hole needle	Box 25
			
AJ 1832A	22g TW x 60mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 21g pre-hole needle	Box 25
			
AJ 1833A	22g TW x 50mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 21g pre-hole needle	Box 25

### 23 GAUGE CANNULA + 21 GAUGE PRE-HOLE NEEDLE

			
AJ 1840A	23g UTW x 70mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 21g pre-hole needle	Box 25
			
AJ 1842A	23g UTW x 50mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 21g pre-hole needle	Box 25

## NEW GTI CANNULA®


### 25 GAUGE GTI CANNULA + 23 GAUGE PRE-HOLE NEEDLE

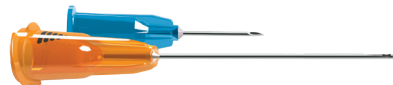			
AJ 1852A	25g UTW x 40mm	GTI Cannula® (Amar) + 23g pre-hole needle	Box 10


Turn to back page for more information about the new GTI Cannula®.

## FAT TRANSFER CANNULAS


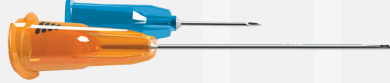
### 25 GAUGE FAT TRANSFER CANNULA

			
AE 1702	18g UTW x 70mm	Straight fat transfer cannula with screw thread hub, closed end and single side port	Box 40


### 20 GAUGE FAT TRANSFER CANNULA

			
AE 1711	20g UTW x 60mm	Straight fat transfer cannula with screw thread hub, closed end and single side port	Box 40


### 25 GAUGE CANNULA + 23 GAUGE PRE-HOLE NEEDLE

			
AJ 1850A	25g UTW x 50mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 23g pre-hole needle	Box 25
			
AJ 1851A	25g UTW x 40mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 23g pre-hole needle	Box 40


### 27 GAUGE CANNULA + 25 GAUGE PRE-HOLE NEEDLE

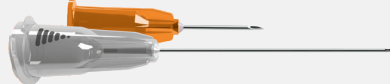			
AJ 1861A	27g UTW x 40mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 25g pre-hole needle	Box 40

### 30 GAUGE CANNULA + 27 GAUGE PRE-HOLE NEEDLE





			
AJ 1870A	30g UTW x 30mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port + 27g pre-hole needle	Box 40

### 18 GAUGE CANNULA

			
AJ 1701	18g UTW x 70mm	Straight dermal cannula with fine thread, two-start hub, closed end and single side port	Box 40

## SHARP NEEDLES


### SHARP NEEDLES

			
M0286	27g TW x 12.7mm	Straight dermal injection needle	Box 100
			
M0288	29g TW x 12.7mm	Straight dermal injection needle	Box 100
			
M0287	30g TW x 12.7mm	Straight dermal injection needle	Box 100
			
M0289A	30g TW x 4mm	Straight dermal injection needle	Box 100

# The latest innovation for small scars and surface defects.

## The new Silkann® GTI Cannula®

Designed specially to disrupt higher density fibrous scar tissue to perform a simple and safe subcision.

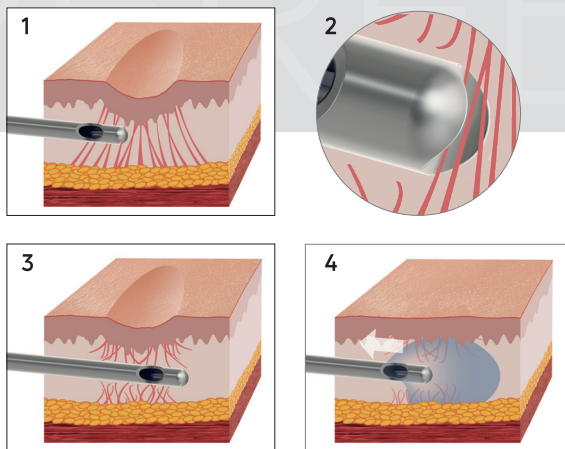

- 1: The Grooved Tipped Injection Cannula is designed to move safely through higher density & fibrotic tissue in the face.
- 2: Micro subcision of fibrous strands.
- 3: Fibrous strands are released.
- 4: Injection of dermal filler to fill the void upon cannula withdrawal.

Designed in conjunction with DR Olivier Amar; an Aesthetic & Reconstructive Surgeon based in London, UK.



**NEW**

**GTI CANNULA®**

**25 GAUGE GTI CANNULA + 23 GAUGE PRE-HOLE NEEDLE**



AJ 1852A	25g UTW x 40mm	GTI Cannula® (Amar) + 23g pre-hole needle	Box 10
----------	----------------	---	-----------

Distributed by:

For more information call:  
+44 (0)1527 501480  
or email: [sales@sterimedix.com](mailto:sales@sterimedix.com)

**Sterimedix®**  
aesthetic

Sterimedix Limited, 1 Madeley Road, North Moons Moat  
Redditch, Worcestershire B98 9NB UK

[sterimedix.com](http://sterimedix.com)