

**pfm**medical

Quality and Experience

---

**Surgical Instruments**

- › Disposable Scalpels
- › Safety Scalpels
- › Surgical Blades and Handles
- › Blade Remover



FEATHER® - World class in quality and consistency for cutting instruments

[www.pfmmedical.com](http://www.pfmmedical.com)




## FEATHER® Surgical Scalpels

### General Benefits

Due to its optimised micrometer-accurate precision grinding, the high-quality stainless steel blade offers outstanding quality and performance. The double-honing process ensures ultra-sharp cutting edges of superb durability.

### General Handling

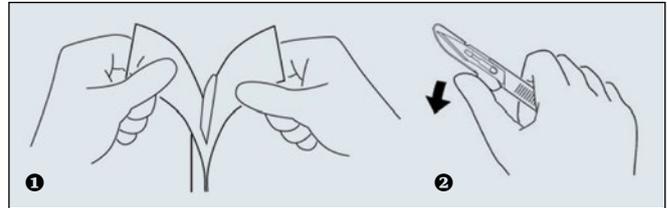

- ❶ **Open the aluminium foil packaging as illustrated:**  
Pull the two flaps apart.
- ❷ **Removal of the blade protection:**  
Press gently on the snap tab as illustrated.

## Standard Scalpels

The disposable Standard Scalpel with high-quality stainless steel blade is suitable for all surgical applications.

### View



### Ordering Information

#### Standard Scalpels

REF	Figuration	PU
200130010	10	20
200130011	11	20
200130012	12	20
205023003	12 D	20
205143003	14	20
200130015	15	20
205353003	15 C	20
200130020	20	20
200130021	21	20
200130022	22	20
200130023	23	20
200130024	24	20
200130025	25	20

## Scalpels, mini

The mini scalpel with its 18 cm handle is ideally suited for all surgical procedures with narrow operating fields.

### View



### Ordering Information

#### Scalpels, mini

REF	Figuration	PU
200170011	11	20
200170014	14	20
200170015	15	20

## Scalpels, long

The scalpel with a plastic handle which is twice the length (23 cm) of a standard model makes it more convenient to use in deep operations.

### View



### Ordering Information

#### Scalpels, long

REF	Figuration	PU
200210010	10	20
200210011	11	20
200210015	15	20
200210020	20	20
200210021	21	20
200210022	22	20
200210023	23	20

## Safety Scalpels


The FEATHER® Safeshield™ disposable scalpel is designed for ultimate safety during surgical procedures.

### View



### Ordering Information

#### Safety Scalpels

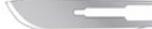REF	Figuration	Pantone	PU
200140010	10 	 471 C	10
200140011	11 	 190 C	10
200140014	14 	 185 C	10
200140015	15 	 265 C	10
200140115	15c 	 2718 C	10
200140020	20 	 465 C	10
200140021	21 	 360 C	10
200140022	22 	 Black 6 C	10
200140023	23 	 3265 C	10
200140024	24 	 405 C	10
200140025	25 	 1585 C	10

### Benefits

#### Optimum safety

The scalpel ensures maximum safety due to its Safeshield™ blade protection which can be closed with one hand.

#### User friendly

Color-coded scalpels for easy blade shape identification

- ▶ Ultra -sharp cutting edges due to high precision grinding technology
- ▶ Single-handed opening and closing of the blade protection
- ▶ Color-coded scalpels for easy blade shape identification
- ▶ Balanced design
- ▶ Outstanding durability
- ▶ Individual sterile packaging

## FEATHER® Surgical Blades and Handles

### General Benefits

The surgical stainless steel blades made from the highest quality hardened stainless steel offer outstanding quality and performance due to their optimised micrometer-accurate precision grinding. The double-honing process ensures ultra-sharp cutting edges of superb durability.

### Easy to use

The blades are ready for use immediately after removal from the packaging.

### Stainless Steel Blades

FEATHER® Stainless Steel Blades made from the highest quality hardened stainless steel are suitable for all surgical applications.

### Ordering Information

#### Stainless Steel Blades

REF	Figuration	PU
201000010	10 	100
201000011	11 	100
201000012	12 	100
205020003	12 D 	100
205140003	14 	100
201000015	15 	100
205350003	15 C 	100
201000020	20 	100
201000021	21 	100
201000022	22 	100
201000023	23 	100
201000024	24 	100
201000025	25 	100

### Carbon Steel Blades

FEATHER® Carbon Steel Blades with an anti-rust coating are suitable for all surgical applications.

### Ordering Information

#### Carbon Steel Blades

REF	Figuration	PU
201500010	10 	100
201500011	11 	100
201500012	12 	100
201500015	15 	100
201500020	20 	100
201500021	21 	100
201500022	22 	100
201500023	23 	100
201500024	24 	100

## Handles

FEATHER® Handles are suitable for carbon steel and stainless steel blades.

### View



### Ordering Information


#### Handles

REF	Version	Figuration	Size	Material	PU
202000003	Standard	10 - 15	3	Stainless steel	10
202000004	Standard	20 - 25	4	Stainless steel	10
202000007	Standard	10 - 15	7	Stainless steel	10
202010003	Long	10 - 15	3L	Stainless steel	10
202010004	Long	20 - 25	4L	Stainless steel	10
208800077	Standard	10 - 15	3DLC	Diamond-like carbon film coating	1
208800078	Standard	20 - 25	4DLC	Diamond-like carbon film coating	1
208800079	Standard	10 - 15	7DLC	Diamond-like carbon film coating	1
208800080	Round	10 - 15	3T	Titanium	1

## Blade Plier

The Feather® Blade Plier is designed to remove sharp surgical blades safely from its handle. It improves safety and hygiene by eliminating direct contact with blades and facilitates replacements.

### View

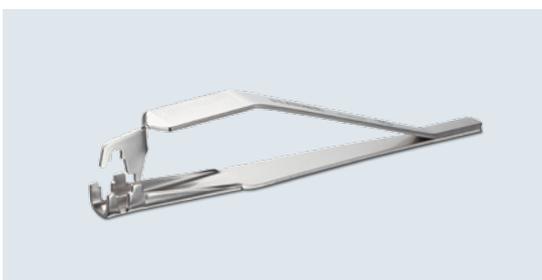

### Benefits

- ▶ Optimal safety
- ▶ Simple „place and pull“ handling
- ▶ Outstanding durability

### Technical Data

- ▶ **For blade figurations:** 10-15
- ▶ **For FEATHER® Surgical Blade Handles:** No. 3, No. 3L, No. 7
- ▶ **Material:** Stainless steel

### Ordering Information

#### Blade Plier

REF	PU
206000002	1

## Blade Remover

### View



Standard Blade Remover



Palm-size Blade Remover

### Benefits

#### Ease of use

The FEATHER® Blade Remover permits the safe and simple removal of used disposable blades from any kind and size of scalpel handle.


#### Safe and hygienic

No direct contact with blades provides the greatest possible protection of medical staff against infection.

#### Disposable product

Once the cost-effective blade remover is filled with used blades, it can simply be disposed of.

### Handling

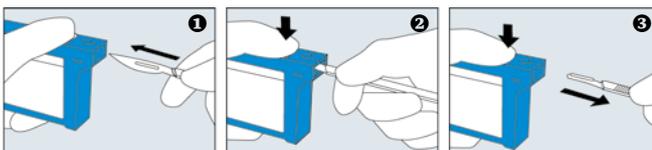

### Instructions for use

- ① Make sure that the correct side of the blade is pointing upwards; then insert the blade into the slot of the blade remover until it hits the stop.
- ② Gently press the button on the top of the blade remover.
- ③ Remove the handle. The removed blade drops safely into the container of the blade remover.

### Ordering Information

#### Blade Remover

REF	Version	PU
201090000	Disposable Blade Remover standard	1
201091000	Disposable Blade Remover palm-size	1

### Contact

Should you have any questions our Customer Service & Sales Support team is happy to advise you.

✉ [service@pfmmedical.com](mailto:service@pfmmedical.com)

☎ +49 (0)2236 9641-220

☎ +49 (0)2236 9641-51

pfm medical ag  
Wankelstraße 60  
50996 Köln, Germany

Certified by  
DIN EN ISO 13485

#### Manufacturer

FEATHER® Safety Razor Co. Ltd., 3-70, Ohyodo-Minami 3-Chome Kita-Ku, Osaka 531-0075, Japan, 0123

Exclusive Distributor for Europe (Except Finland, Switzerland)

pfm medical ag, Wankelstr. 60, 50996 Köln, Germany

[www.pfmmedical.com](http://www.pfmmedical.com) - Your information source for the products of pfm medical ag.

PB2304EN/2018-11-27

